


# Ambu® Intubation Trainer

45

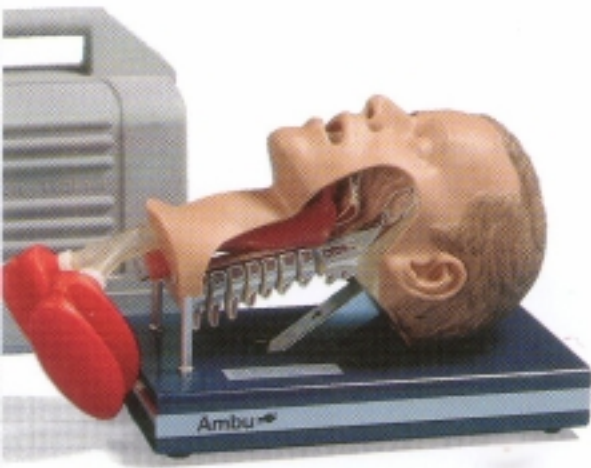

The Ambu Intubation Trainer is for teaching intubation techniques with laryngoscopes, airways, endotracheal tubes, LMA's, Combitube and other auxiliary aids for airway management.

The accurate simulation of the mouth, nostrils, teeth, tongue, pharynx, epiglottis, vocal cords, trachea, oesophagus and lungs, coupled with realistic movement of the head, cervical spine and jaw simulating relevant anatomical changes during intubation, makes it the most realistic Intubation Trainer on the market.

The left side of the head is cut away so a student's performance can be supervised and the walls of the pharynx and trachea are transparent enabling the student to follow the tube down the throat.

Acoustic signals triggered by excess pressure on the teeth allows the student to learn when mistakes are made and how and when to correct them. If desired the sensitivity of this signal can be adjusted.

Different patient types can be practiced on by adjusting the difficulty of intubation.

## TECHNICAL DATA

Size: 40 x 24 x 35 cm (with transport cover) :  
Following tubes are recommended:

Weight: 9 kg.

Endotracheal tube 8 mm, Nasotracheal tube 7 mm,  
Laryngeal mask and LMA-fastrach, sizes 3 and 4

## ORDERING INFORMATION


## CAT. NO.

Ambu Intubation Trainer in hard transport case

000 186 000

# Ambu® IV Trainer

45

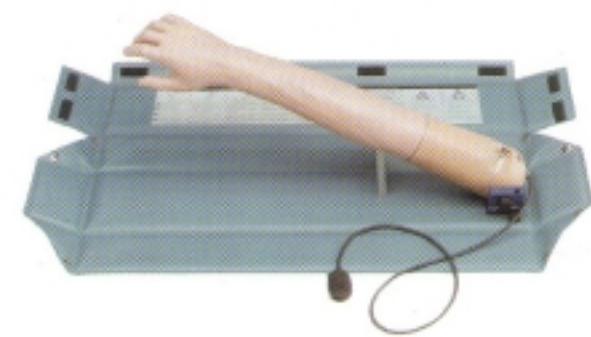

The Ambu I.V Trainer is for training the insertion of cannulas and catheters, the infusion of fluids and the injection of medication, blood sampling and pulse measurement.

## Ambu IV Trainer Features

- Replaceable skin and veins are made of natural latex, which features hole closure after penetration, extending the lifetime of the arm
- The movable wrist and rotating arm ensures an absolutely realistic training product, especially when a catheter is placed on the back of the hand

- The valve system reduces the formation of air pockets in blood vessels
- The indicator panel which, represents each blood vessel, gives immediate feedback e.g. if an infusion has been correctly accomplished
- The multi-purpose carry case serves both as a washable underlay and as a protective carry case. It also includes a bracket so that an infusion stand can easily be mounted



## TECHNICAL DATA

Size (L): of arm 68 cm  
Carry case dimensions: 75 x 15 x 11 cm

Reservoir capacity: 260 ml.

Weight: 3.0 kg

The I.V. Trainer can be mounted on all Ambu Man models



## ORDERING INFORMATION

## CAT. NO.

Ambu IV Trainer, complete with infusion stand, red dye concentrate, talcum powder and 5 ml and 50 ml syringes.  
Supplied in a carry case/washable underlay

255 001 000